

Medtrica Mattresses Models




	Standard	Symmetry	Symmetry Plus	Memory Flex 3.0	Memory Flex 5.0
Construction	Single Core	Dual Layer	Dual Layer; Two Zone - Body and Heel	Triple Layer, 4 Zone - Head/ Shoulders, Torso, Calf, Heel	Triple Layer, 4 Zone - Head/ Shoulders, Torso, Calf, Heel
Firm Internal Side Rails	No	Yes	Yes	Yes	Yes
Memory Foam Layer	No	No	No	Yes - 3.0 lb Viscoelastic	Yes - 5.0 lb Viscoelastic
User Weight Rating	300 lbs	450 lbs	450 lbs	550 lbs	550 lbs
Warranty	5 year Non Prorated on Cover and Core	5 year Non Prorated on Cover and Core	5 year Non Prorated on Cover and Core	7 year Non Prorated on Cover and Core	7 year Non Prorated on Cover and Core
Fluid Proof Cover Fabric options	70D Nylon - (no butyl backing); 4-way Stretch Dartex™	70D Nylon - (no butyl backing); 4-way Stretch Dartex™	70D Nylon - (no butyl backing); 4-way Stretch Dartex™	70D Nylon - (no butyl backing); 4-way Stretch Dartex™	70D Nylon - (no butyl backing); 4-way Stretch Dartex™
Anti-Shear Inner Liner	Optional	Optional	Optional	Optional	Optional
Low Slip Bottom	Yes	Yes	Yes	Yes	Yes
4 Reinforced Handles	Yes	Yes	Yes	Yes	Yes
Fire Rating	All materials used exceed flammability standards Cal TB 117, 16CFR 1632	All materials used exceed flammability standards Cal TB 117, 16CFR 1632	All materials used exceed flammability standards Cal TB 117, 16CFR 1632	All materials used exceed flammability standards Cal TB 117, 16CFR 1632	All materials used exceed flammability standards Cal TB 117, 16CFR 1632
Anti-Bacterial, anti-Fungal, anti-Static, Stain & Tear resistant. Latex Free	Yes	Yes	Yes	Yes	Yes
100% Medical Grade Foam and Fabric	Yes	Yes	Yes	Yes	Yes
Interface Pressure Score Results	NA	50% Pressure Relief; 50% Pressure Reducing	75% Pressure Relief; 25% Pressure Reducing	75% Pressure Relief; 25% Pressure Reducing	75% Pressure Relief; 25% Pressure Reducing
Braden Scale	4	3 and 4	2, 3 and 4	2, 3 and 4	2, 3 and 4
No Cost Options	No Zippers or Handles	Sloped Heel Section; No Zippers or Handles; No Internal Side Rails	Sloped Heel Section; No Zippers or Handles; No Internal Side Rails	Sloped Heel Section; No Zippers or Handles; No Internal Side Rails	Sloped Heel Section; No Zippers or Handles; No Internal Side Rails
Options	Extra Fire Barrier to meet Cal TB 129 & BFC IX-11; Sloped Heel Section, Anti-Shear Liner, Firm Internal Side Rails; Custom Sizes	Extra Fire Barrier to meet Cal TB 129 & BFC IX-11; VersaCare™ Style; Custom Sizes	Extra Fire Barrier to meet Cal TB 129 & BFC IX-11; VersaCare™ Style; Custom Sizes	Extra Fire Barrier to meet Cal TB 129 & BFC IX-11; VersaCare™ Style; Custom Sizes	Extra Fire Barrier to meet Cal TB 129 & BFC IX-11; VersaCare™ Style; Custom Sizes
Standard Sizes	76" x 35" x 6" 78" x 35" x 6" 80" x 35" x 6" 84" x 35" x 6"	76" x 35" x 6" 78" x 35" x 6" 80" x 35" x 6" 84" x 35" x 6"	76" x 35" x 6" 78" x 35" x 6" 80" x 35" x 6" 84" x 35" x 6"	76" x 35" x 6" 78" x 35" x 6" 80" x 35" x 6" 84" x 35" x 6"	76" x 35" x 6" 78" x 35" x 6" 80" x 35" x 6" 84" x 35" x 6"
Bariatric Weight Rating and Warranty	NA	950 lb user Weight Rating, 5 year non-prorated cover on core and cover	950 lb user Weight Rating, 5 year non-prorated cover on core and cover	1000 lb user Weight Rating, 5 year non-prorated cover on core and cover	1000 lb user Weight Rating, 5 year non-prorated cover on core and cover

Medtrica Mattresses

Medtrica offers a full range of 100% medical grade foam and fabric mattresses.

From the Single Layer Standard Mattress to the 3 Layer, 4 Zone SolaCare Memory Flex 5.0 model, Medtrica has a mattress that is the **perfect No Risk fit** for your requirements regarding:

- Patient Safety
- Patient Comfort
- Facility Budget

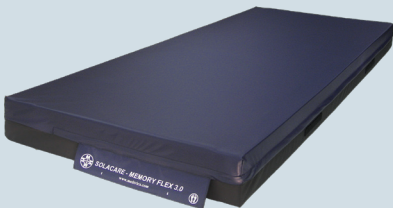

Medtrica mattresses provide optimum value because they are engineered for **maximum patient safety and comfort while easily fitting into the common budget concerns of your facility.**

HIGHER QUALITY FOAM MEANS:

- Weight ratings that **support a much greater weight range of user** without the added cost of specialty bariatric models.
- Lower interface pressure scores for **greater patient safety and lower facility risk.**
- **Mattresses hold their shape and performance capabilities** throughout their life-time.

Firm internal side rails directly address FDA entrapment guidelines.

70D Nylon fabric without a butyl backing **provides superior patient comfort and safety** because it creates less heat build-up between the user and the surface. It also provides a better fit and is far more durable / lasts longer than covers made from 70D nylon fabric with butyl backing.

Dartex™ 4-way stretch covers are the highest rated in the industry because they **provide excellent interface pressure scores** and also **do not breakdown** when cleaned with standard or harsh disinfectants like other polyurethane 4-way stretch covers do.

Anti-Shear liners create less friction and heat build up between the user and the mattress thus creating a more comfortable, cooler and safer surface.

Exceed OEM bed frame manufacturer's mattress specifications.

For Acute and Long Term Care

30 Day Trial Available
Full Non-prorated warranties on Cover and Core
Custom Sizing and Construction Available

Manufactured by:

MEDTRICA SOLUTIONS LTD.
6350 204th Street, #104A
Langley, BC. Canada V2Y 2V1
t: 604-618-2291
f: 604-952-4314



info@medtrica.com
www.medtrica.com

MEDTRICA SOLUTIONS LTD

4001 Inglewood Ave
101-128
Redondo Beach, CA 90278
t: 877-703-1079
f: 310-263-9098